

Agilent BioTek 406 FX
Microplate Washer Dispenser
User Manual



Notices

Document Identification

2121000N

Revision C. TBD 2024

Copyright

© Agilent Technologies, Inc. 2021-2023

No part of this manual may be reproduced in any form or by any means (including electronic storage and retrieval or translation into a foreign language) without prior agreement and written consent from Agilent Technologies, Inc. as governed by United States and international copyright laws.

Agilent Technologies, Inc.
5301 Stevens Creek Blvd.
Santa Clara, CA 95051
USA

Software Revision

This document is valid for the 3.0 revision or higher of the Agilent Liquid Handling Control (LHC) and compatible programs until superseded.

Instrument Manufacturing



Manufactured by
Agilent Technologies, Inc.
5301 Stevens Creek Blvd.
Santa Clara, CA 95051
USA



Warranty

The material contained in this document is provided "as is," and is subject to being changed, without notice, in future editions. Further, to the maximum extent permitted by applicable law, Agilent disclaims all warranties, either express or implied, with regard to this manual and any information contained herein, including but not limited to the implied warranties of merchantability and fitness for a particular purpose. Agilent shall not be liable for errors or for incidental or consequential damages in connection with the furnishing, use, or performance of this document or of any information contained herein. Should Agilent and the user have a separate written agreement with warranty terms covering the material in this document that conflict with these terms, the warranty terms in the separate agreement shall control.

Technology Licenses

The hardware and/or software described in this document are furnished under a license and may be used or copied only in accordance with the terms of such license.

Restricted Rights Legend

U.S. Government Restricted Rights. Software and technical data rights granted to the federal government include only those rights customarily provided to end user customers.

Agilent provides this customary commercial license in Software and technical data pursuant to FAR 12.211 (Technical Data) and 12.212 (Computer Software) and, for the Department of Defense, DFARS 252.227-7015 (Technical Data - Commercial Items) and DFARS 227.7202-3 (Rights in Commercial Computer Software or Computer Software Documentation).

Safety Notices

CAUTION

A **CAUTION** notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in damage to the product or loss of important data. Do not proceed beyond a **CAUTION** notice until the indicated conditions are fully understood and met.

WARNING

A **WARNING** notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in personal injury or death. Do not proceed beyond a **WARNING** notice until the indicated conditions are fully understood and met.

Contents

Notices	i
1 General Information	1
Contact Information	2
Worldwide Sales and Support	2
Technical Support and Service	2
Customer Care	2
Intended Use Statement	2
Quality Control	2
Safety Notices	2
Warnings and Precautions	3
Electrical Hazards	3
Chemical/Environmental	4
Components	5
Intended Product Use	6
Safety Symbols	6
Conformance to Standards and Directives	12
Directives	12
Standards	12
EMC Information and Technical Description	13
Emissions	13
Immunity	13
Disposal	13
2 Introduction	14
Introducing the 406 FX Microplate Washer Dispenser	15
Package Contents	16
Accessory Pack	16
Optional Accessories	16
General Instrument Accessories	16
Washer Accessories	17
Washer Manifolds	18
Magnetic Bead Assay Accessories	18
Miscellaneous	18
Peri-pump Optional Accessories	19
Peri-pump Modules	19
Cassettes and Accessories	19
Tips and Cassette Tubing Extension Kits	20
General accessories	20

Syringe Dispenser Optional Accessories	21
Manifolds	21
Bottles and Tubing	21
General accessories	21
Integration Kits	22
Features of the 406 FX	22
Washer	22
Peri-pump Dispenser	23
Syringe Dispensers	23
Liquid Handling Control (LHC) Software	23
BioStack Compatibility	24
406 FX Dispenser Comparison	25
3 Installation	27
Important Information	28
Unpack the Instrument	28
Remove the Shipping Hardware	29
Supplies	29
Setting Up the 406 FX	30
Connect the Vacuum Pump, Tubes, and Bottles	31
Waste System	31
Install the Vacuum Line Filter	33
Final Check	33
Attach the Mist Shield	33
Install the Microplate Carrier	34
Install the Left-Side Carrier Kit (optional)	35
Set Up the Peri-pump Dispensers	36
Install the Secondary Peri-Pump	36
Dispense Cassette Diagram	39
Install the Dispense Cassette	40
Change the cassette type setting	42
Prime Trough Inserts	42
Install the Syringe Dispenser Component	44
Install Inline Filter for 32-Tube Dispensers	45
Install Tubing and Manifolds for Syringe Dispenser	46
Check the Syringe Manifold Installation	50
Syringe Dispenser Check Valves	50
Install Software/Connect to Computer	50
Connect to Host Computer	51
Connect Power Cable	51

Define Instrument Settings	51
LHC Users Only	51
Standard Vacuum Pump Users	52
Define Startup Preferences (LHC users)	52
Verify Performance	52
Verify the Peri-Pump Dispenser	53
Verify the Syringe Dispensers	53
Verify the Washer Component	54
Repackaging and Shipping Instructions	55
Repacking the Syringe Dispenser	61
4 Operation	63
Basic Operation	64
Optimize Performance	64
Keep the devices clean and the tubing wet	64
Prime the tubing to remove air bubbles	64
Best Practices for all 406 FX Devices	65
Before using the Peri-pump	65
Before using the Syringe Dispenser	66
Optimize protocols to improve evacuation	66
Optimize protocols to specify the desired residual volume	66
RUN: Running Predefined Protocols	67
Creating Protocols: Washing, Aspirating and Dispensing Fluid	67
Touch Screen Control	67
Liquid Handling Control (LHC) Software	67
Customize the Predefined Protocols	68
Predefined Protocols	69
Predefined Sample Protocols	70
Protocol Parameters Tables	72
Syringe Dispenser Flow Rates	74
Wash Step Parameters Table	75
Operating with the BioStack	77
Touch Screen Control:	77
LHC Control:	78
Touch Screen Basics	79
Quick Prime	80
Quick Prime	80
Quick Wash	80
Quick Wash	80
Learn more:	80

Quick Clean	80
Quick Clean (Sonicate)	80
Add protocols to the Home screen	81
Remove protocols from the Home screen	81
Create, Edit, or Copy a Protocol	82
Copy a Protocol	83
Lock (password-protect) and unlock protocols	84
Rename a protocol	84
Change the Plate Type	85
Using LHC to Control the 406 FX Washer Dispenser	85
Predefined Protocols	85
Communications Port	85
406 FX Steps	86
LHC Protocols	86
Customize the Protocols	86
LHC Protocols Explained	87
Prerequisite	87
Protocol Files	87
Managing Multiple Instruments	88
Validate versus Run	88
Target Instrument Settings	89
Get settings from:	89
How to define a Protocol (LHC only)	90
Customize the Predefined Protocols	91
How to define the Plate Map	91
Plate Types and Processing Patterns	92
Plate Types Table	93
Handling special plates	93
1536F - 1536-well Flanged Plates	93
Handling Special Plates	94
1536F - 1536-well Flanged Plates	94
About Wash Processing Patterns	95
Dispense Processing Patterns	97
Columns	97
Rows	97
Peri-pump Dispense Pattern	98
Dispense Pattern with Both Syringes	102
Washer Operation	102
About the 406 FX Wash Step	102

Wash Step Parameters Table	102
Shake/Soak Step Parameters	105
Pre-Dispense	105
Washer Flow Rates	106
Recommendations for 1536-well hardware	107
How to define a Final Aspiration step	108
Add a secondary aspiration to a wash cycle	108
Cell Wash	109
Low-flow Fluid Path	109
Define a Cell Wash Protocol	110
Cell Wash Strategies	111
Best practices:	112
Biomagnetic Separation - Magnetic Bead Assays	113
Handling and Cleaning the Magnets	113
Realign the BioStack with the Magnet Installed	114
Perform Magnetic Bead Assays	114
Optimize Magnetic Bead Protocols	114
Flat Magnets	115
Ring Magnets	117
How to determine the Magnet Adapter Height Offset	118
Magnet Height Offset	118
Bottom-Wash	119
Washer Flow Rates	120
Peri-pump Peristaltic Dispenser	121
Quick Dispense-Prime-Purge	121
Dispense Cassettes	122
Peri-pump: How it works	122
Recommended prime volumes for the Peri-pump	123
At the start of the day:	123
Release the tension on the dispense cassette	124
Peri-pump Dispense Step	125
Define the Dispense Step	125
Require a Specific Peri-pump Cassette	125
Shorten the Dispense Cassette Tubing	126
Standard Cassette Dead Volume	127
Peri-pump Advanced Settings	127
Settings for 0.5 μ L Dispense	127
Dual Syringe Dispenser	128
Quick Dispense	128

How to Prime the Syringe dispenser	128
Recommended prime volumes for the Syringe dispensers	130
Syringe Prime Step	130
Syringe Dispense Step	131
Change the Syringe Dispenser Manifold	132
Special Procedure for Magnetic Bead Assays	133
406 FX Instrument Utilities: General	135
Adjust Utility	137
Positional Limits	137
How to use the Adjust Utility- LHC	137
Use the Align Tool to optimize manifold tube positioning	138
AutoPrime	139
About AutoPrime	139
Specify the Interval and AutoPrime Parameters	139
Turn on AutoPrime for the Washer	140
Turn on AutoPrime for the Peri-pump Dispenser	140
Turn on AutoPrime for the Syringe Dispenser	141
Delay After Dispense	142
Standard Vacuum Pump Users	142
Upload-Download Protocols (LHC Only)	143
Operating with the BioStack	144
Touch Screen Control:	144
LHC Control:	144
BioStack Alignment Utility	146
Plate Lid Definition	146
Transfer Plate Lid Definitions	147
Use the BioStack	147
Re-stack plates in the BioStack	148
Nth Plate to perform special process	149
Using Nth Plate to perform special functions	149
About the Nth Plate Step	151
Advanced Plate Sequencing with the BioStack	152
5 Maintenance	157
Preventative Maintenance Overview	158
Washer	158
Peri-pump	158
Syringe Dispenser	158
Recommended Maintenance Schedule	159
Recommended Maintenance Schedule	160

Daily Maintenance	161
Overnight/Multi-Day Maintenance	162
AutoPrime for Overnight-Weekend Maintenance	162
Soak the manifold tubes in cleaning fluid	164
Prime Step	164
Create a Prime protocol	165
Removing Protein Residuals and Fungi Growth	166
Daily Practice with buffer and deionized water:	166
Weekly or As Needed, use NaOH and HCl to remove proteins:	166
Alternatively, use an Enzyme-Active Detergent:	167
Direct Drain Only: Clean the Waste Bottle Filter	167
Periodic Maintenance	167
Clean the Bottles	168
Clean the Plate Carrier	169
Clean the Exterior Surfaces and Mist shield	169
Autoclavable Components	169
Washer Maintenance	170
AutoClean the Washer	170
Create an AutoClean protocol	171
Touch screen	172
AutoClean the 32-Tube Dispense Manifold	172
Tools for cleaning dispense and aspirate tubes	173
Remove and clean the washer manifold	173
Clean the Buffer Bottle Filter	175
Peri-pump Dispenser Maintenance	175
Flush the Dispense Cassette	176
Tools and Supplies	176
Unclog the Dispense Tips	177
Record the Number of Plates Processed	178
Cassette Expected Lifetime	178
Syringe Dispenser Maintenance	178
Clean the Bottles and Tubes	179
Syringe Dispenser Maintenance Tasks	179
Clean the Manifolds	179
Clean the Dispense Tubes	180
Clean or Replace the Check Valves	181
Autoclave the Syringe Head	182
Decontamination	183
Tools and Supplies	183

Decontaminate Exterior Surfaces	184
Decontaminate Tubing and Manifold	184
Alternate Decontamination Procedure for Tubing and Manifold	187
Long Shutdown (Prepare for Storage or Shipment)	188
Storing the Instrument	189
Replace Components	189
Replace Washer Manifold O-rings and Seals	189
Replace Dispense Cassette Tips	191
Calibrate the Backlash for Syringe Dispenser	192
Plate-Type-Specific Syringe Manifolds	192
6 Qualification	194
Qualification Overview	195
Washer	195
Washer Qualification Tests	195
Peri-pump and Syringe Dispensers	195
Qualification Schedule	196
Liquid Testing the 406 FX	196
Which Tests to Perform?	196
Washer	196
Peri-pump Dispenser	197
Syringe Dispenser	197
Important Recommendations for All Liquid Tests	197
System Self-Test, Verify Information	198
Record Onboard Software Information	198
Sensor Test: Fluid, Waste	198
Test the Fluid Supply Sensor	198
Test the Waste Level Sensor	199
Ultrasonic Cleaner Test	199
Washer Dispense Precision and Accuracy Test	199
Predefined Test Protocols	200
Predefined Protocols	200
Washer Qualification Test Materials	201
Evacuation Efficiency Test (96)	202
Evacuation Diagnostic Test	203
Dispense Precision and Accuracy Test (96)	204
Buffer Switching Dispense Test	205
Evacuation Efficiency Test (192)	205
Evacuation Diagnostic Test (192)	206
Dispense Precision and Accuracy Test (192)	207

Evacuation Efficiency Test (128)	208
Dispense Precision and Accuracy Test Worksheet	209
96-Tube Manifold	209
Evacuation Efficiency Test Worksheet	210
96-Tube Manifold	210
Dispense Precision Test Worksheet 96-Tube Manifold Annual Buffer Switching Test	211
Dispense Precision Test Worksheet 192-Tube Manifold	212
Evacuation Efficiency Test Worksheet 192-Tube Manifold	213
Dispense Precision Test Worksheet 192-Tube Manifold Annual Buffer Switching Test	214
Evacuation Efficiency Test Worksheet 128-Tube Aspirate Manifold	215
Peri-pump Dispense Precision and Accuracy Tests	216
Dispense Precision and Accuracy Specifications	216
Peri-pump Dispenser Test Materials	217
Peri-pump Dispenser Precision and Accuracy Test Solutions	217
Perform the Peri-Pump Precision and Accuracy Tests	218
Syringe Dispenser Liquid Tests	220
Syringe Dispenser Test Materials	220
Syringe Dispenser Test Solutions	220
Perform the Syringe Dispense Precision & Accuracy Test	221
Syringe Dispenser Precision & Accuracy Test Worksheet Test 1 / 96-Well Microplate / 40 μ L Dispense	223
Syringe Dispenser Precision & Accuracy Test Worksheet Test 1 / 384-Well Microplate / 20 μ L Dispense	224
Syringe Dispenser Precision & Accuracy Test Worksheet Test 2 / 96-Well Microplate / 160 μ L Dispense	225
Syringe Dispenser Precision & Accuracy Test Worksheet Test 2 / 384-Well Microplate / 80 μ L Dispense	226
Syringe Dispenser Precision & Accuracy Test Worksheet Test 3 / 1536-Well Microplate / 6 μ L Dispense	227
Syringe Dispenser Precision & Accuracy Test Worksheet Test 4 / 96-Well Microplate / 300 μ L Dispense	228
7 Troubleshooting	229
Troubleshooting	230
Error recovery	230
Error Codes	230
General Operation Troubleshooting	230
Startup	231
Microplate Scratches	231
Cell-based Assays	231

Washer Troubleshooting	231
Fluid Aspiration	232
Fluid Delivery	233
Fluid Leakage	235
Washer Manifold Movement	235
Vacuum Pump Problems	236
Peri-pump Troubleshooting	236
Syringe Dispenser Troubleshooting	237
Fluid Delivery	237
Fluid Leakage	239
Dispense Manifold Movement	239
Microplate Carrier Movement	239
Syringe Movement	240
Communication Errors	240
Appendix A: Peri-pump Test Worksheets	242
Peri-pump Dispense Precision & Accuracy Test Worksheet1 μ L Test for 1 μ L Cas- sette	243
Peri-pump Dispense Precision & Accuracy Test Worksheet10 μ L Test for 1 μ L Cas- sette	244
Peri-pump Dispense Precision & Accuracy Test Worksheet1536 Test	245
Peri-pump Dispense Precision & Accuracy Test Worksheet5 μ L Test	246
Peri-pump Dispense Precision & Accuracy Test Worksheet50 μ L Test for 5 μ L Cas- sette	247
Peri-pump Dispense Precision & Accuracy Test Worksheet 10 μ L Test for 10 μ L Cas- sette	248
Peri-pump Dispense Precision & Accuracy Test Worksheet100 μ L Test	249
Appendix B: Specifications	250
Physical Specifications	251
Performance Specifications	253
Washer	253
Peri-Pump Dispenser	253
Cassette Expected Lifetime	254
Syringe Dispensers	254
Appendix C: Error Codes and Conditions	256
Error Codes Overview	257
When an Error Code or Error Message Appears:	257
System Error Codes	258
406 FX-LHC Software Error Codes	266

Appendix D: Changing the Washer Manifold	272
Changing the Wash Manifold	273
To reconfigure the instrument, first update the washer manifold setting:	274
Storing the Unused Manifolds	274
Appendix E: Safety Information	275
Safety Notices	276
Warnings and Precautions	278
Electrical Hazards	278
Chemical/Environmental	280
Components	287
Intended Product Use	290
Chemical Compatibility Chart	294
In This Book	298

DRAFT



2 Introduction

This chapter introduces the 406 FX and provides instructions for contacting Technical Support.

Introducing the 406 FX Microplate Washer Dispenser	15
Package Contents	16
Optional Accessories	16
Peri-pump Optional Accessories	19
Features of the 406 FX	22
BioStack Compatibility	24
406 FX Dispenser Comparison	25

DRAFT

2 Introduction

Introducing the 406 FX Microplate Washer Dispenser

Introducing the 406 FX Microplate Washer Dispenser

The 406 FX offers up to three devices in one instrument: a microplate Washer, a peristaltic pump dispenser (the "Peri-pump"), and a dual Syringe pump dispenser. One or both dispensers are provided with every 406 FX.

One model is available: 406FX--SN.

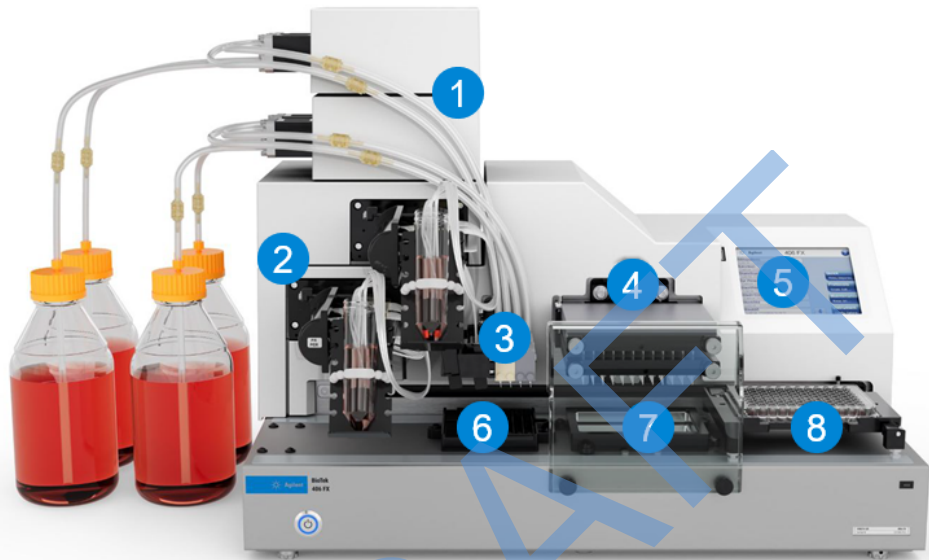


Figure 2-1: 406 FX Washer Dispenser

Device/Component	Description
1 Syringe Dispensers	Syringe-pump dispensers, each with an 8-, 16- or 32-tube manifold. Up to two syringe pumps may be connected.
2 Peri-pump Dispensers	Peristaltic, 8-channel dispensers with an entirely visible fluid path. Shown with optional secondary Peri-pump.
3 Dispense arm	Holds the Peri-pump Tip Holder and the Syringe dispenser manifolds.
4 Washer manifold	96- or 192-tube devices aspirate and dispense fluid; or 128-tube manifold aspirates fluid from 1536-well plates.
5 Touchscreen	On-board instrument control.
6 Priming trough	Waste reservoir for collecting priming fluid for the Peri-pump and Syringe dispensers.
7 Ultrasonic cleaning and priming trough Mist shield	All models support AutoClean ultrasonic cleaning of the washer.
8 Plate Carrier	Holds standard microplates for processing.

Package Contents

Description	Part number
Power supply, 24 V DC	02285
External Vacuum Pump Switch Box (with control cable)	2120505
USB cable	75108
Microplate carrier	2120508
Mist shield and 2 thumbscrews	2122074
Reagent bottle holder	7212034
Reagent bottle holder step, 2 count.	7212035
Phillips screwdriver for removing dispenser arm shipping bracket	01188
Syringe priming trough	2122043
Peri-pump priming trough, 2	7182043
406 FX User Manual, on USB flash drive	2121000N

Accessory Pack

Description	Part number
Hex wrench, 9/64" for removing washer manifold shipping bracket	48434
Hex wrench, 1/16" for removing magnets from syringe manifolds and for adjusting manifold stop screw	48713
Phillips screwdriver 0.3 mm OD, for Peri-pump cassette calibration	98268
Disposable syringe, 10 cc Luer, for Peri-pump cassette maintenance	49919
Stylus, .012 OD no color, for cleaning 192-tube dispense manifold	7772029
Stylus, .018 OD yellow, for cleaning washer and Syringe manifold dispense tubes, and 128-tube aspirate manifold tubes	7772030
Stylus, .031 red, for cleaning washer manifold aspirate tubes	7772033
Manifold stop height set tool for adjusting manifold stop screw	2122136
Shipping set	2123008

Optional Accessories

General Instrument Accessories

Description	PN
Agilent liquid testing solutions for instrument qualification tests	Wetting Agent 7773002 Blue Test Dye 7773001
Liquid Handling Control Software 3.0	5330209

2 Introduction

Optional Accessories

Washer Accessories

Description	Part Number
Complete dispense/waste system 4 L bottles	2120531
Complete dispense/waste system High Flow 4 L bottles	2120532
Complete dispense/waste system High Flow 10 L bottles	2120533
Complete dispense/waste system High Flow 20 L bottles	2120534
Vacuum System Complete Assembly 115V/230V, 4 L bottles	7100746S
High Flow Vacuum System Complete Assembly 115V/230V, 4 L bottles	7100753S
Complete Direct Drain Waste System, 115V/230V, 4L bottles	1170560S
Complete Direct Drain High Flow Waste System 115V/230V, 4 L bottles	1170562
Direct drain waste system upgrade kit, 115V/230V, 4 L bottles (For use with any complete waste system)	1170561S
Complete dispense system. 4 L, Assembly	2120535
Complete dispense system. 10 L, Assembly	2120536
Complete Waste System 115V/230V, 4 L bottles	7100746S
High Flow Vacuum System, Complete Assembly 115V/230V, 4 L bottles	7100753S
Vacuum Set with Level Sense Assembly, 4 L	7100543S
Vacuum Set with Level Sense Assembly, 10 L	7100582S
Vacuum Pump Assembly 115V/230V	7103035
High Flow Vacuum Pump Assembly 115V/230V	7100754

Repackaging and Shipping Instructions

Please read the information provided below before preparing the instrument for shipment.

- Contact Technical Support before returning equipment for service.
- Decontamination prior to shipment is required by the U.S. Department of Transportation regulations.
- If the instrument has been exposed to potentially hazardous material, decontaminate it to minimize the risk to all who come in contact with the instrument during shipping, handling, and servicing. The Maintenance section contains decontamination instructions.
- Install the shipping hardware (see next section).
- The instrument's packaging design is subject to change. If the instructions in this document do not apply to the packaging materials you are using, contact Technical Support for guidance.
- Be sure to use packaging materials supplied by the manufacturer. Other forms of commercially available packaging are not recommended and can void the warranty.
- If the packaging materials have been damaged or lost, or if the same set has been used more than four times, order replacements.

CAUTION

Shipping Hardware. All shipping hardware must be removed before operating the instrument and reinstalled before repackaging the instrument for shipment.

Prepare the instrument for shipping:

LHC	Touch screen
<ol style="list-style-type: none">1 Select Tools > Instrument Utilities > Washer2 Select the shipping bracket link under Service Functions.	<ol style="list-style-type: none">1 Select Instrument > Next > Next > Other.2 Select Park Manifold.
<ol style="list-style-type: none">3 Remove the Peri-pump's cassette and the Syringe dispensers' manifolds.4 Wrap rubber bands around the Peri-pump pump heads5 Remove the plate carrier and put it in the accessories box.6 Uninstall and repack separately the Syringe pump and its accessories, if applicable.7 Slide the transport rail into position next to the dispense arm. It will be secured in place by one of the shipping brackets.	

Reinstall the shipping hardware by reversing the steps described to [Remove the Shipping Hardware](#).


NOTE

The instrument's packaging design is subject to change over time. If the instructions in this section do not appear to apply to the packaging materials you are using, please contact Agilent Technical Support for guidance.

3 Installation

Repackaging and Shipping Instructions


- 1 Repack the instrument accessories box bottom tray, placing each item into its designated slot.

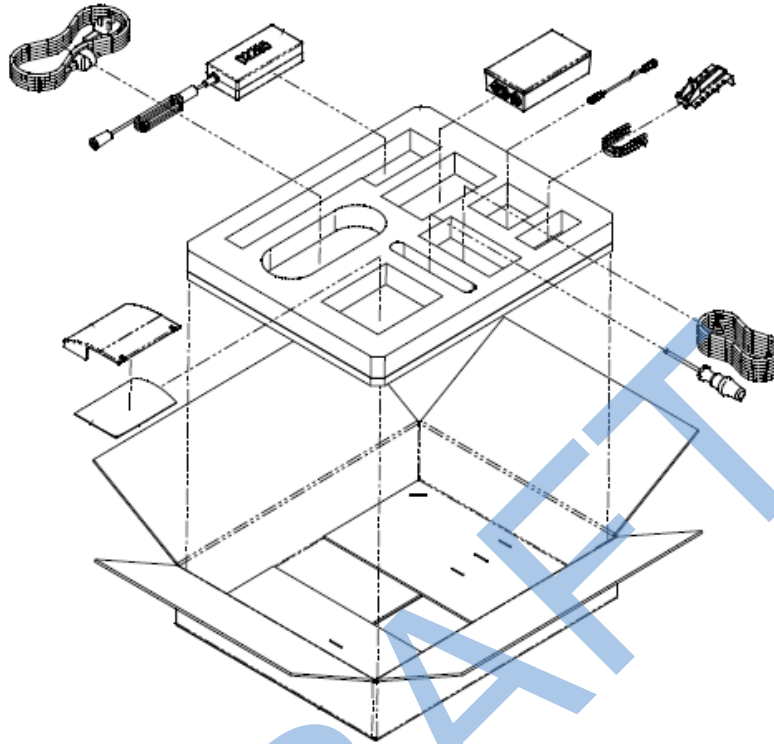

3 Installation

Repackaging and Shipping Instructions

- 2 Repack the instrument accessories box top tray, placing item into its designated slot.
- 3 Place this tray in the box and cover with the tray cover.

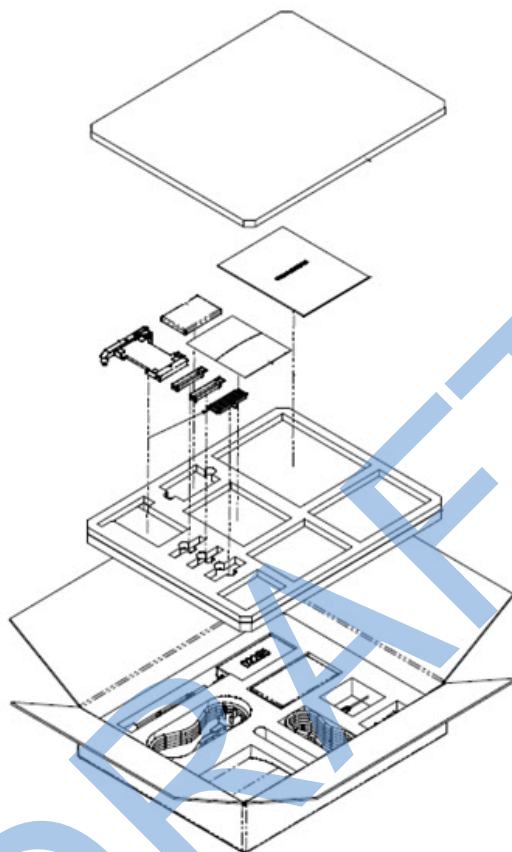


- 4 Close the flaps and tape shut.
- 5 Prepare the main box for the instrument and accessories

3 Installation

Repackaging and Shipping Instructions


- 6 Place the instrument box (A) on a flat, sturdy surface and insert the bottom (B).
- 7 Insert the box insert (C) so that the handholds of the insert align with the outer box. This insert will not fill the entire space.
- 8 Place the prepared accessories box (D) between the insert and the box as see below.

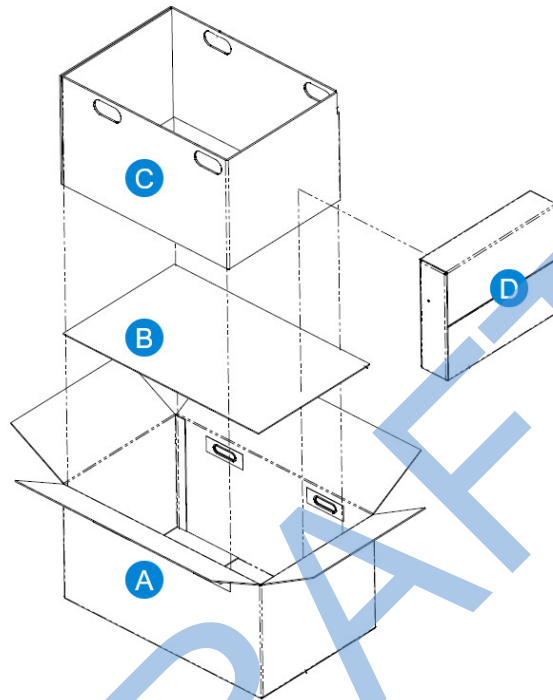

3 Installation

Repackaging and Shipping Instructions

- 9 Place the bottom shipping foam insert into the box and then place the bottom into the foam insert. These support the instrument during shipping.

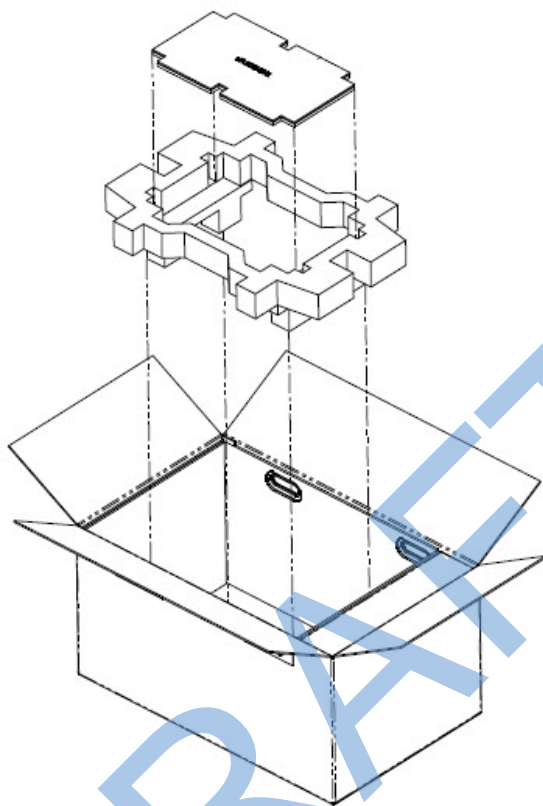


IMPORTANT

The instrument is heavy. Have another person help you with the next packing steps to avoid damage to the instrument.

3 Installation

Repackaging and Shipping Instructions

- 10 Place the instrument in the poly bag.
- 11 Carefully place the instrument in the box, fitting it into the bottom foam insert.
- 12 Place the top foam insert on top of the instrument and then place the cardboard over it.

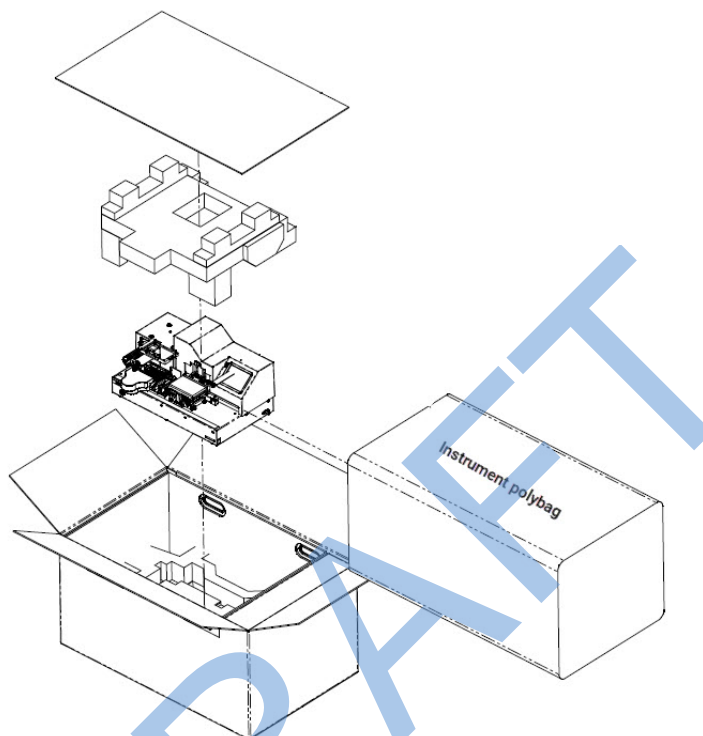


- 13 Securely tape up the box for shipment.

3 Installation


Repacking the Syringe Dispenser

Repacking the Syringe Dispenser

After preparing the instrument for shipping, and reversing the installation steps, pack the Syringe dispenser as shown here:

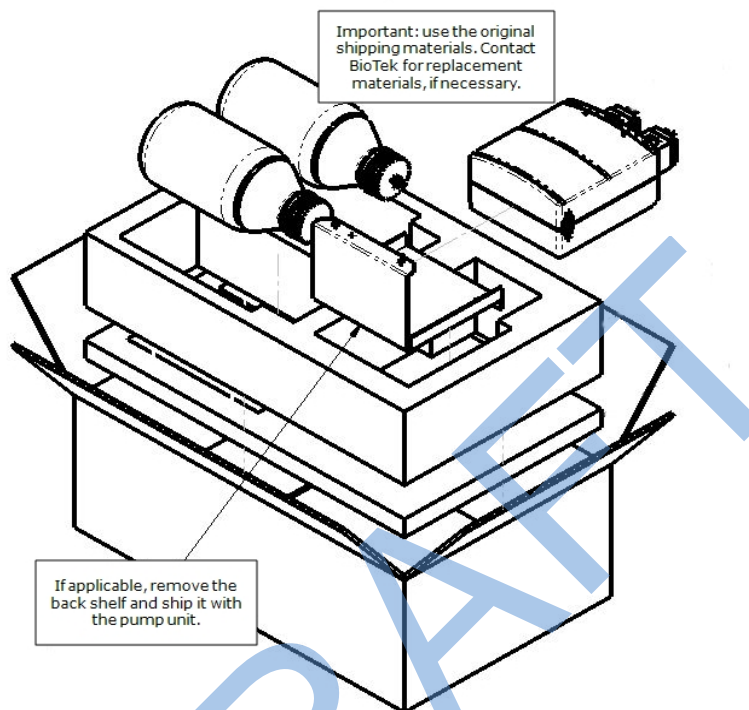

Fill the foam tray with the Syringe dispenser accessories and place the tray on top of the inner foam box containing the pump and supply bottles.

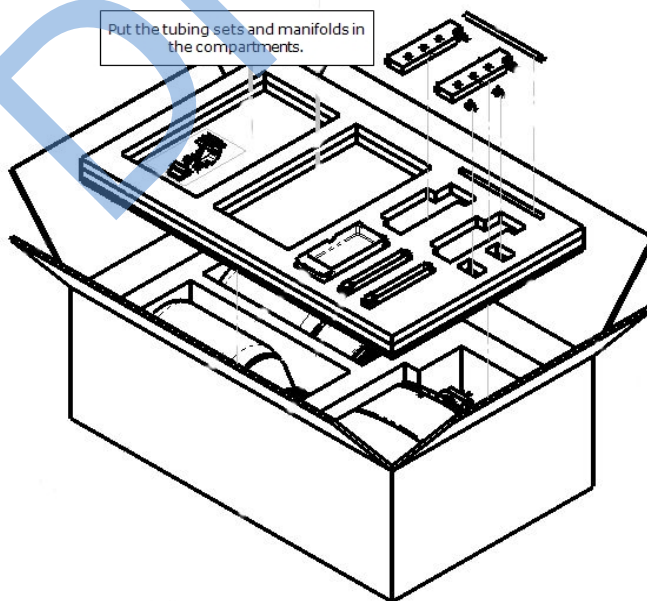

5 Maintenance

Properly maintaining the 406 FX is the key to reliable performance.

Preventative Maintenance Overview	158
Recommended Maintenance Schedule	160
Daily Maintenance	161
Overnight/Multi-Day Maintenance	162
Soak the manifold tubes in cleaning fluid	164
Removing Protein Residuals and Fungi Growth	166
Periodic Maintenance	167
Washer Maintenance	170
AutoClean the Washer	170
Remove and clean the washer manifold	173
Peri-pump Dispenser Maintenance	175
Syringe Dispenser Maintenance	178
Decontamination	183
Long Shutdown (Prepare for Storage or Shipment)	188
Replace Components	189
Calibrate the Backlash for Syringe Dispenser	192

Removing Protein Residuals and Fungi Growth

CAUTION

Bovine Serum Albumin. Solutions containing proteins, such as bovine serum albumin (BSA), will compromise instrument performance over time unless a strict maintenance protocol is adhered to.

When using protein solutions or similar fluids, Agilent recommends performing the following procedures to thoroughly flush out protein particles and other contaminants from the fluid path.

Also note some components can be autoclaved to sterilize them.

NOTE

Four-liter volumes specified in the following are approximate amounts.

Daily Practice with buffer and deionized water:

If the 406 FX is idle between plates for longer than 45 minutes, flush the proteins using buffer. Use buffer for the AutoPrime, as well, for immediate reuse of the instrument. Otherwise, if the tubes will be allowed to dry, flush the system with DI water to remove the buffer:

- 1 Fill a supply bottle with buffer solution. Connect the bottle to the Buffer A port or to the syringe dispenser.
- 2 Run the applicable **DAY_RINSE** protocol.
- 3 Enable AutoPrime for 60-minute intervals.

NOTE

Flush the buffer from the system using DI water when the device will be idle for an extended period, that is, repeat steps 1 and 2 using water.

At the end of the day:

- 1 Fill a supply bottle with buffer for each device (washer or syringe). Connect the bottle to the Buffer A port for the washer or to the appropriate Syringe dispenser port.
- 2 Run the applicable **DAY_RINSE** protocol three times.
- 3 Fill a supply bottle with deionized water and repeat steps 1 and 2.
- 4 Perform your regular Overnight/Multi-Day Maintenance routine.

Weekly or As Needed, use NaOH and HCl to remove proteins:

- 1 Flush the system with 0.1-0.5 N* NaOH (sodium hydroxide), followed by neutralization with an equivalent normality (0.1-0.5 N) of HCl (hydrochloride).
- 2 Rinse well with deionized water to remove the HCl.
- 3 Run the applicable **W-DAY_RINSE** protocol three times with deionized water if you plan to use the device immediately.

NOTE

* N = Normal solution, which contains 1 'gram equivalent weight' (gEW) of solute per liter of solution. The gram equivalent weight is equal to the molecular weight expressed as grams divided by the 'valency' of the solute.

5 Maintenance

Periodic Maintenance

Alternatively, use an Enzyme-Active Detergent:

- 1 Mix an enzyme-active detergent according to the manufacturer's directions to fill two four-liter supply bottles.
- 2 Run the **W-DECONTAMINATE** or **S-DECONTAMINATE** protocol, as appropriate, based on the syringe channel.

NOTE

Follow the on-screen prompts and messages. **Delay** messages contain instructions and information to help you walk through the procedure.

- 3 Respond to the **Delay** message and when ready, press **Continue**.
- 4 When the protocol is completed, determine if the appropriate Decontaminate protocol needs to be run again and repeat as necessary.
- 5 When the process is completed, connect a bottle containing four liters of deionized water and run **W-DAY_RINSE** or **S-DAY_RINSE** as appropriate, based on the syringe channel, three times to flush the system.
- 6 Repeat the procedure for the other Syringe dispensers as appropriate.


Direct Drain Only: Clean the Waste Bottle Filter

Periodically clean the direct-drain waste bottle filter (part number 01310):

WARNING

Direct Drain Waste Bottle Filter: Handle with care! When cleaning the waste bottle filter, be aware of any potential biohazards and handle appropriately.

- 1 Open the bottle and lift the cap and its tubing up and out of the bottle.
- 2 Remove the fluid filter at the bottom of the tubing.
- 3 Wash the filter with hot water and a soft-bristle brush, if necessary.
- 4 Rinse the filter and reinstall it.

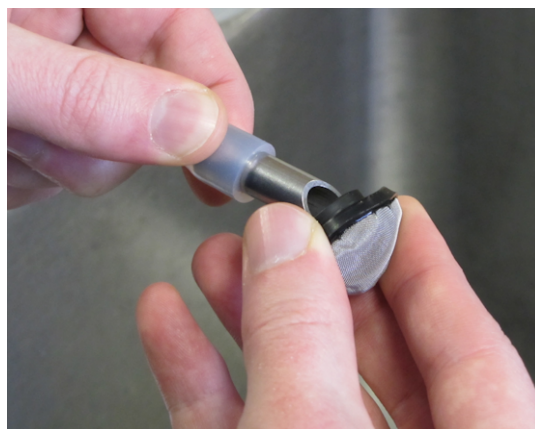

Periodic Maintenance

Periodic maintenance involves cleaning the components on a regular basis to keep the instrument running efficiently and in compliance with performance specifications. The recommended **frequency for cleaning components** is *at least monthly*. The risk and performance factors associated with your assays may require that some or all procedures be performed more frequently.

WARNING

Internal Voltage. Always turn off the power switch and unplug the power supply before cleaning

5 Maintenance

Periodic Maintenance

the outer surface of the instrument.

WARNING



Potential Biohazards. Some assays or specimens may pose a biohazard. Adequate safety precautions should be taken as outlined in the assay's package insert. Always wear safety glasses and appropriate protective equipment, such as chemical-resistant rubber gloves and apron.

WARNING

Liquids. Avoid spilling liquids on the instrument; fluid seepage into internal components creates a potential for shock hazard or instrument damage. If a spill occurs while a program is running, stop the program and turn off the instrument. Wipe up all spills immediately. Do not operate the instrument if internal components have been exposed to fluid.

CAUTION

Liquids. Do not immerse the instrument, spray it with liquid, or use a dripping-wet cloth on it. Do not allow water or other cleaning solution to run into the interior of the instrument. If this happens, contact Technical Support. Do not soak the touch screen.

CAUTION

Lubricants. Do not apply lubricants to moving parts. Lubricant on components in the carrier compartment will attract dust and other particles, which may cause the instrument to produce an error. Do not apply lubricants to manifold o-rings, channel-end seals, bottle cover seals, any tubing connection, or any surface that is a part of the fluid path.

CAUTION

Chemical Compatibility. Some chemicals may cause irreparable damage to the instrument. The following chemicals have been deemed safe for use in the instrument: buffer solutions (such as PBS), saline, surfactants, deionized water, 70% ethyl, isopropyl, or methyl alcohol, and 20% sodium hydroxide. Never use acetic acid, DMSO, or other organic solvents. These chemicals may cause severe damage to the instrument.

CAUTION

Caution: Bovine Serum Albumin. Solutions containing proteins, such as bovine serum albumin (BSA), will compromise instrument performance over time unless a strict maintenance protocol is adhered to.

Perform these maintenance tasks regularly:

- [Clean the Bottles below.](#)
- [Clean the Plate Carrier on the facing page.](#)
- [Peri-pump Dispenser Maintenance on page 175.](#)
- [Syringe Dispenser Maintenance on page 178.](#)

Clean the Bottles

- Clean and rinse the supply bottles with deionized water before the first use, before each refill, and as necessary, to prevent bacteria growth.
- Empty the waste bottle often (at least daily), and firmly seat the waste bottle stopper.
- Rinse the covers every time the wash or rinse bottles are filled.
- Accumulated algae, fungi, or mold may require decontamination.
- Check the hex nuts securing the quick-disconnects to the bottle cap to ensure they are not loose or corroded.

NOTE

To ensure that fluid does not back up into the vacuum pump during operation, always operate the instrument with the **waste sensor cable** installed and the **waste detection sensor** enabled (the

5 Maintenance

Periodic Maintenance

sensor is enabled by default).

NOTE

If fluid collects in the **overflow bottle**, thoroughly rinse the level-switch assembly and bottle.

Clean the Plate Carrier

If liquid has overflowed onto the plate carrier, transport rail, or glide strips, some buildup may occur and prevent the microplate from seating correctly on the carrier. This buildup can interfere with plate transport. Weekly cleaning is recommended.

- 1 Turn the instrument off.
- 2 Lift the carrier up and off the transport rail.
- 3 Clean the carrier, rails, and glide strips using mild detergent and hot water, 70% isopropyl alcohol, or ethanol. Clean the priming trough, too.
- 4 If detergent was used, wipe the components with a cloth moistened with water. Use a clean, dry cloth to dry the components.
- 5 Reinstall the carrier:
 - Place it on the transport rail so the slot on its base fits into the **Y-axis Carrier Arm**.
 - If necessary, release the spring-loaded microplate clamp in the back left corner of the carrier to level the carrier on the base.

Clean the Exterior Surfaces and Mist shield

- 1 Turn off the instrument and disconnect the power cable.
- 2 Moisten a lint-free disposable towel with water or with water and mild detergent. **Do not soak the cloth.**
- 3 Remove the mist shield if it is attached. Wipe the inside and outside surfaces of the mist shield with the towel. Wipe the top surface of the instrument base and all exposed surfaces of the instrument.
- 4 If detergent was used, wipe all surfaces with a water-moistened cloth.
- 5 Use a clean, dry cloth to dry all wet surfaces.

Autoclavable Components

Autoclaving is an efficient method of sterilizing instrument components. It is a good alternative to some of the decontamination procedures for qualified items.

Do autoclave:	Do NOT autoclave:
Peri-pump cassette*	Washer manifolds
16-Tube and Low-density plate-type-specific Syringe dispenser manifolds	Plate carrier
Priming trough inserts	32-Tube Syringe dispenser manifolds
Syringe pump head (glass/stainless steel)	8-Tube Syringe dispenser manifolds (gray block holds two manifolds, 16 tubes total)

Chemical Compatibility

For your own safety and to ensure the instrument's long life, follow these guidelines when choosing compounds to use.

Important: Avoid continuous contact with harsh chemicals. A best practice is to rinse all components with deionized/distilled water after contact with any strong acid, base, or solvent.

Materials listed as 1-20 in **Table 1** are cross-referenced as headings 1-20 in **Table 2**. Table 1 indicates which component contains the material using **P** for Peri-pump, **S** for the dual Syringe dispensers, and **W** for the washer module.

Table 1: Material Where Used

#	Material		Where Used
1	304 Stainless Steel	W	Inlet screen, feeder tube to manifold, vacuum switch, fluid pump, direct-drain inlet screen
2	316 Stainless Steel	P	Steel dispense tips, pump rollers, carrier rails, cassette bolts
		S	Dispense tubes, syringe piston, cylinder head and fittings, check valve spring, aspirate tubes, aspirate arm release button
		W	Fluid inlet filter, feeder tube to manifold, vacuum switch, fluid pump, and vacuum filtration pressure regulator components
3	Aluminum		Instrument main structural plate
		P	Center Holder, Tube Tensioner, pump head, dispense arm
4	Borosilicate Glass	S	Syringe cylinder, supply bottles (autoclavable models only)
5	Brass	P	Threaded inserts in Center Holder and Tip Holder
6	EP (Ethylene Propylene)	S	Check-valve o-rings (organic solvent tubing set)
		W	Inlet valve, vacuum switch
7	ETFE (Ethylene tetrafluoroethylene)	S	Manifold fitting
8	Neoprene	W	Manifold channel-end seals
9	Nylon	W	Inlet fitting, vacuum switch adjustment screw, carrier leveling feet
9	Polyethylene	S W	Supply Bottle
10	Polycarbonate		Plate carrier
		P	Plate positioning dog, peristaltic pump cover latch
		W	Vacuum switch
11	Polyethylene	W	Buffer bottle
12	Polypropylene	P	Molded dispense tips
		S	Bottle cap and fittings, check valves (organic solvent tubingset)
		W	Outlet fitting, fittings in bottles, inline fittings, float ball, vacuum filtration pressure gauge components, direct-drain module fittings
13	Polystyrene		Assay plates, mist shield
		P	Peristaltic pump cover, RFID antenna cover
		W	Flow sensor
14	PTFE (polytetrafluoroethylene)	S	Syringe seal
		W	Optional check valves (PN: 68098) for fluid pump
15	CPVC (Polyvinyl chloride)	W	Inlet screen, inlet valve, manifold, waste sensor, ball retainer, vacuum filtration carrier, system orifice, and pressure regulator, direct drain bottle sensor

#	Material	Where Used	
16	Ryton PPS	P	Tube Organizer, priming trough, priming trough insert, RAD 8-to-1 chute
		S	Manifold, manifold plugs, priming trough insert
		W	Outlet valve, fluid pump, inlet valve, direct-drain fluid pump
17	Santoprene	W	Fluid pump, direct-drain pump check valves
18	Silicone	P	Tubing
		S	Tubing, manifold o-rings, check valve o-rings
		W	Inlet tubing, outlet tubing, o-rings, vacuum filtration pressure regulator tubing and o-rings, direct-drain module tubing
19	Ultem (polyetherimide)	P	Center Holder, Tube Tensioner, Tip Holder
		S	Check-valves
		W	Outlet valve, cell wash inlet valve, vacuum filtration valve, direct-drain valve
20	Viton	W	Outlet valve, direct-drain valve

DRAFT

[illegible]

Chemical	1-20 represent Materials listed above in the Where Used List																			
	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20
Sodium Hypochlorite <20%	C	D	A	D	B	A	A	C	A	A	A	A	A	C	B	B	A	C	B	B
Sodium Hypochlorite 0.50%	B	D	A	D	B	A	A	C	A	A	A	A	A	C	B	A	A	C	B	A
Sulfuric Acid <10%	B	B	A	N	A	C	A	A	A	A	A	A	A	A	C	A	A	A	C	A
Trichloroethylene	B	D	N	N	N	A	D	N	D	C	D	A	D	A	D	D	D	A	D	D
Virkon 10%	A	D	A	N	A	N	A	A	A	A	A	N	A	A	A	N	A	A	A	N

* Acetone is not recommended for use with Agilent instruments.

† Exposure to DMSO and Acetonitrile may cause the silicone tubing to swell, increasing the volume of fluid Dispensed. The magnitude of this effect will vary with concentration and exposure. Re-calibration of the cassette may be required.